HOMEDICS®

Register your product today at
www.homedicsgroup.com/register

Microdermabrasion

Instruction Manual

MDA100-EU

2 year
guarantee

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

IMPORTANT SAFEGUARDS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- This appliance can be used by children aged from 14 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre or similarly qualified person for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- This appliance is for household and indoor use only.
- This product requires a 220-240 V AC power.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, send to a HoMedics Service Centre. All servicing of this appliance must be performed by authorised HoMedics service personnel or similarly qualified person only.
- Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch. Keep air openings free of lint, hair etc.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the mains supply when removing the battery. The battery is to be disposed of safely. The consumer must use the appropriate tools to remove the screws from the product, then disengage the clip fittings on the outer casing to access the battery for removal, prior to the product being scrapped.

SAFETY CAUTIONS:

PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.


- **If you have any concerns regarding your health, consult a doctor before using this appliance.**
- If you are susceptible to allergic reactions or if you have circulatory problems, consult a physician before using this product.
- Do not use on irritated, sunburned or chapped skin, or on any area with dermatological problems.
- Do not apply over skin eruptions, warts, moles or varicose veins.
- Use should be easy and painless. If you experience skin irritation, or hand discomfort when using this product, discontinue use. If a more extreme reaction occurs, consult your doctor
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by Diabetics
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- If the product feels excessively hot, switch off and contact the HoMedics Service Centre.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP.
- Ensure that your hands are dry when operating the controls or removing the plug.
- This is a non professional appliance and is not intended to treat and medical conditions of the skin. Do not use as a substitute for medical attention.
- This product should NEVER be used by any individual suffering from any ailments that would limit the user's capacity to operate the controls.

Microdermabrasion is NOT recommended if you have any of the following conditions:

- Lesions, warts, or recent herpes outbreak
- Sunburn or sun sores
- Active rosacea
- Auto-immune disorders
- Skin cancer
- Vascular lesions
- Any open wounds, sores, or cuts
- Taking oral blood thinners (anti-coagulants)

PRODUCT FEATURES:

HoMedics Microdermabrasion brings professional salon skincare into your own home. The diamond tip works with a gentle but powerful suction action to remove dead cells from the skin surface, revealing and stimulating the fresh skin layer beneath. Perfect for restoring youthful radiance to the skin; improving skin tone and texture, smoothing the appearance of fine lines and wrinkles and reducing the appearance of age spots and sun damage. Use regularly to maintain young and fresh looking skin.

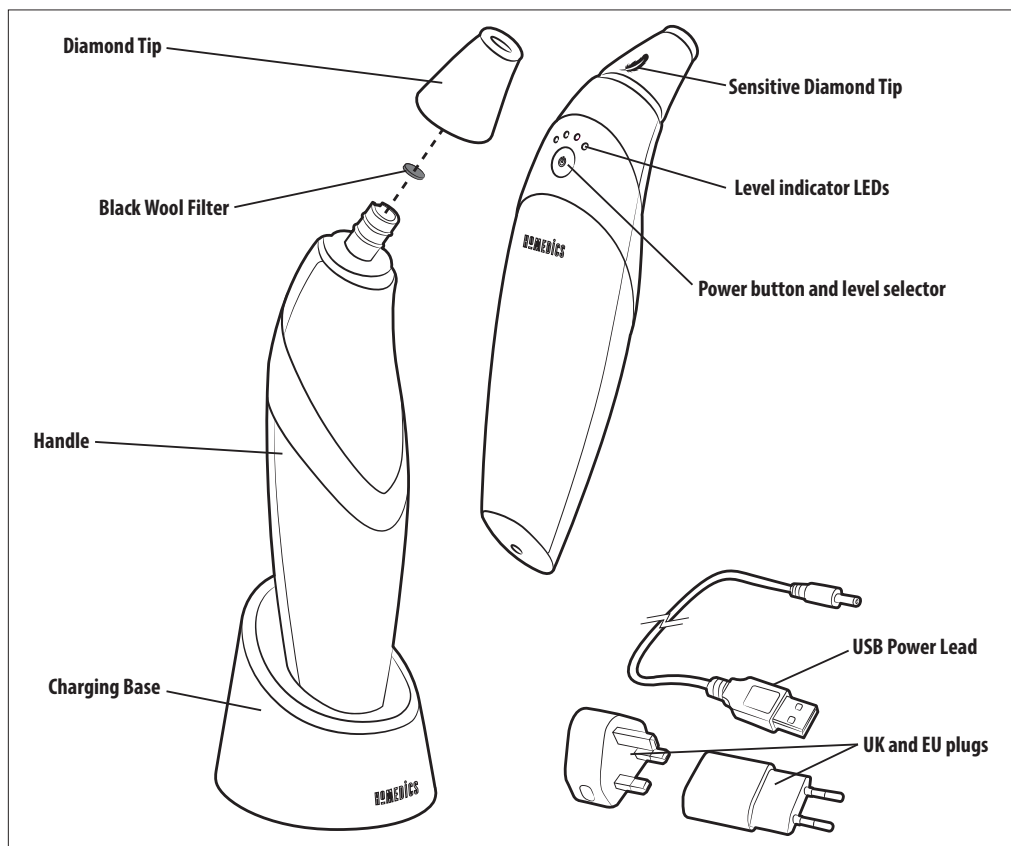Battery directive

 This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

MAINTENANCE

To Store

Place attachments into their proper compartments inside the storage case. Store case in a safe, dry place, out of children's reach.

To Clean

Always unplug before cleaning the unit. Never use liquids or abrasive cleaners to clean the housing of the unit or the case. Wipe off case with a damp cloth. Clean attachments periodically using detergent. Allow attachments to dry before inserting them into the hand unit.

Adapters


UK Model Number: SSA021F050100BSUM

EU Model Number: SSA021F050100EUUM

Input: 100-240V ~50/60Hz 0.2A

Output: 5V  1A

INSTRUCTIONS FOR USE:

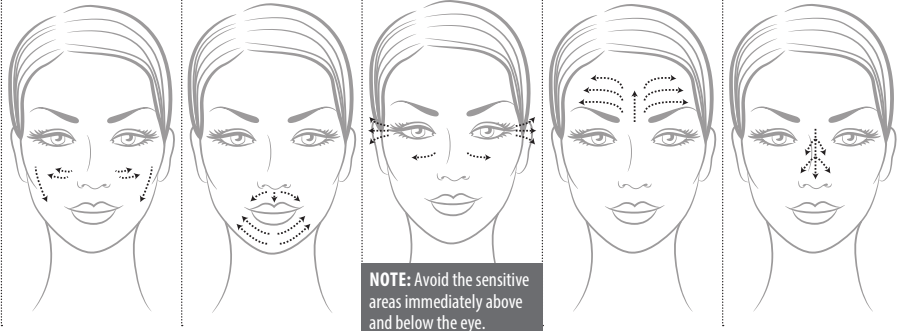AREA	Cheeks	Around Lips & Chin	Around Eyes	Forehead	Nose
MAX LEVEL	3	2	1	4	4
ROUTINE	 NOTE: Avoid the sensitive areas immediately above and below the eye.				

Charging the battery

When the level indicator LED flashes red, the product needs charging. Connect the charging base to the mains using the USB power lead and the appropriate plug. Place the handle in the charging base and press down firmly. The LED lights will flash to indicate that charging is taking place. When charging is complete all the lights will stay on. It takes 4 hours for a full charge, which will provide 90-120 minutes use, dependent on power levels used.

Filter Fitting and Replacement

For correct operation it is very important that a new black wool filter is fitted before each use of the device. The filter collects dead skin cells and should be disposed of after use and then replaced. Using the device without a filter will allow debris and dead skin cells to enter the motor, shortening the life of the product and invalidating the guarantee. (Fig.1)

What to Expect


Regular microdermabrasion treatments can both restore and maintain a youthful appearance to the skin. We advise starting slowly with just one or two treatments a week, then building up to treatments every other day if desired.

During use you should feel a gentle suction effect as the tip moves across the skin and some people may also notice the slight coarseness of the diamond tip. Depending on your skin type, you may experience some reddening or increased sensitivity after treatment. This is perfectly normal and usually subsides within a few hours. You may wish to consider this and make your first treatment at the end of the day rather than just before going out.

Adverse reactions are rare, but if skin becomes unduly inflamed or irritated, stop using the device immediately.

Facial Treatment Routine

- Wash and dry your face
- Remove the handle from the charging base and fit your choice of diamond tip with a new black wool filter.
- Press and hold the power button for 2 seconds to turn the unit on. Press again as required to increase the power level, as indicated by the LEDs. It is recommended to start on level 1 for the first treatment and to progress to a higher level if comfortable. See also MAX recommended levels in the above diagrams

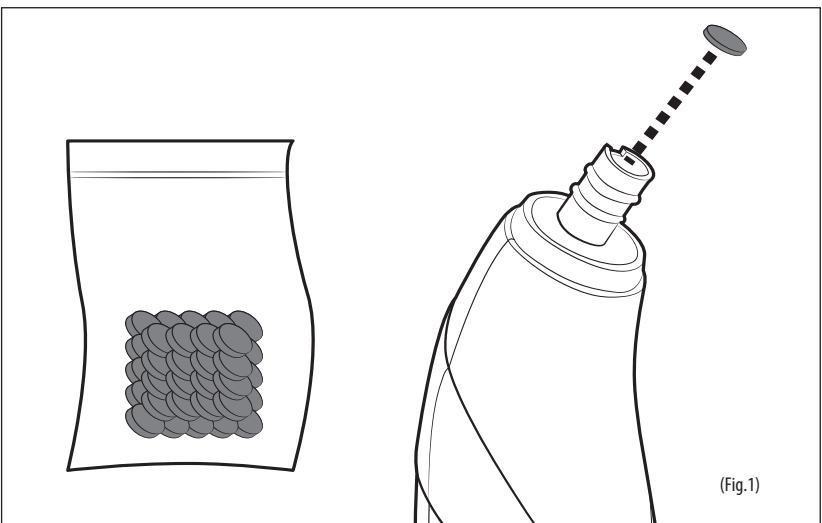(Fig.1)

- Place the diamond tip flat against the skin until you feel suction then move it slowly across the surface of the skin in a smooth motion, covering the same area no more than 2-3 times in one treatment. Do not press hard. In some areas it can help to hold the skin taught with the other hand (eg. around the eyes).
- The diagrams above will help you establish an effective facial routine.

After treatment

- Press and hold the power button for 2-3 seconds to turn off device.
- Wash your face with cold water to remove any remaining dead skin cells and then apply your preferred moisturiser. **NOTE:** Avoid using AHA (acid based) moisturisers on the day of treatment
- Depending on your skin type, you may experience some reddening or increased sensitivity after treatment. This is perfectly normal and usually subsides within a few hours.
- Avoid direct sun exposure after treatment, consider applying a stronger sunscreen as required.
- Always, dispose of the used black wool filter and replace with a new filter ready for next use.
- Remove dead skin cells and debris from the diamond tip using the brush provided. We also recommend the periodic use of an ultrasonic cleaner for deep cleaning the diamond tip, especially if it starts to feel smooth or less effective.

Other Treatments

Other areas of the body can also be treated 2-3 times a week. For example, neck, chest, hands etc. When treating larger areas, replace the black wool filter part way through the treatment in order to maintain the level of suction. Aim to cover the same area no more than 2-3 times in one treatment.

FAQ

How often should I use the product?

It's best to start slowly – usually just once a week, then build up to 2 or 3 times a week if this is comfortable for your skin.

Will it hurt?

No. Treatment should not be painful. Depending on the sensitivity of your skin, you may feel slight discomfort from the diamond tip as it moves across the skin and there could be some short term redness or sensitivity after treatment, but this will soon subside. If treatment is uncomfortable, then you may need to use a lower power setting or switch to the sensitive diamond tip.

Can I use the product in combination with anti-aging creams/serums?

If your anti-aging cream or serum contains AHA (or acid) then it should not be used on the same day as the microdermabrasion treatment, but can be used on non-treatment days. Non-AHA (acid free) anti-aging products can be used at any time.

The product does not switch on

The power button has to be pressed continuously for 2-3 seconds to turn the product on. If you try this and it still doesn't work, then it may need recharging. See 'CHARGING THE BATTERY'

The product seems to be losing power

This can happen if the black wool filter is not replaced after each treatment. Try fitting a new filter. Also check the tip to ensure that the plastic part of the tip has no cracks or damage. Check that the two rubber o-rings are in place on the nozzle end of the unit (where the tip fits over).

The diamond tip no longer seems rough

The diamond tip can accumulate skin cells over time which may start to make it feel smooth and reduce its effectiveness. If this happens we recommend cleaning the diamond tip in an Ultrasonic Cleaner. These are either available to purchase, or jewellery stores offer a cleaning service. If using the product regularly you may need to do this several times a year to maintain optimum effectiveness.

What if I still have a problem or question that is not answered here?

Please contact our service team who will be very happy to help you.

Customer Service

Tel: +44 (0) 1732 378 557

Email: support@homedics.co.uk

www.homedics.co.uk

ACCESSORIES AND SPARE PARTS

Available from our website: www.homedics.co.uk/parts-and-accessories.html

- Replacement Black wool filters (50pcs)
- Diamond tip
- Sensitive diamond tip
- Precision tip - for precise treatment of small areas
- Body tip - for use over larger areas.

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